

Specification Sheet

Liquid Chromatograph Mass Spectrometer

LCMS-8050



Triple Quadrupole Mass Spectrometry is the method of choice for definitive identification and reproducible quantification of trace-level analytes in complex samples for a variety of applications including clinical research, forensic toxicology, pharmacokinetics, environmental analysis, and food and beverage testing. Combined with the chromatographic resolving power of our world-leading UHPLC systems, and maintaining Shimadzu's proprietary ultrafast technologies (UF Technologies), which include high-speed MRM transitions, MS/MS acquisition, and ultra-high speed polarity switching, the LCMS-8050 can dramatically improve analytical throughput with ultra-high-speed performance. In addition, the newly designed ion source and collision cell, Heated-ESI and UFsweeper® III collision cell technology, yield higher sensitivity and expand the potential range of LC/MS/MS applications.

Instrument

Model	LCMS-8050
Mass range	<i>m/z</i> 10 to 2000
Sensitivity	ESI positive: 1 pg reserpine, S/N > 60,000:1 (RMS)
Resolution	R < 0.7 u FWHM
Mass stability	0.05 u/12hr
Cross-talk	< 0.003 %
Minimum pause time	1 msec
Minimum dwell time	0.8 msec
Scan speed	Max 30,000 u/sec (in all modes of scanning) (0.1 u step: 300,000 data points/sec)
Polarity switching time	5 msec
Interface	ESI (Standard), APCI (Optional), DUIS (Optional)
Applicable LC flow rate	ESI 1 µL/min to 2 mL/min
MRM transition speed	Max 555 channels/sec
DL maximum temperature	300 °C
Block heater maximum temperature	ESI/DUIS: 500 °C APCI: 300 °C
Interface maximum temperature	ESI/DUIS: 400 °C APCI: 500 °C

Analysis mode	Q1 Scan/SIM Q3 Scan/SIM MRM Precursor ion scan Product ion scan Neutral loss scan
---------------	--


Mass Analyzers and Detector

Mass analyzers	Q1 & Q3 are molybdenum hyperbolic mass filters with pre-rods; Q1 includes post-rods
Collision cell	Tapered multipole type ultra-high-speed collision cell (UFsweeper® III collision cell)
Detector	Secondary electron multiplier with off-axis conversion dynode
Ion optics	Q-array focus optics operating in Field-Flow mode, multipole transfer optics
Digital detection system	Operates in pulse counting mode for fastest operation
Detection mode	Ultra-fast positive / negative ion switching
Vacuum system	Rotary pump 1 unit Vacuum pumping speed 28 m ³ /hour Triple-inlet turbo molecular pump 1 unit 40 L/sec, 260 L/sec, 210 L/sec

Software

Workstation	LabSolutions LCMS Version 5.6 for LCMS-8050
Operating system	Windows 7 32 Bit, 64 Bit
Instrument control	Prominence and Nexera series LCMS-8050 and interface
MS acquisition mode	Scan (Max. 1,000 events), SIM (Max. 1,000 events × 32 channels)
MS/MS acquisition mode	MRM (Max. 1,000 events × 32 channels) Product ion scan Precursor ion scan Neutral loss scan
Auto-tuning	Possible to optimize sensitivity, resolution, and mass calibration in both positive and negative ionization mode

Installation Example

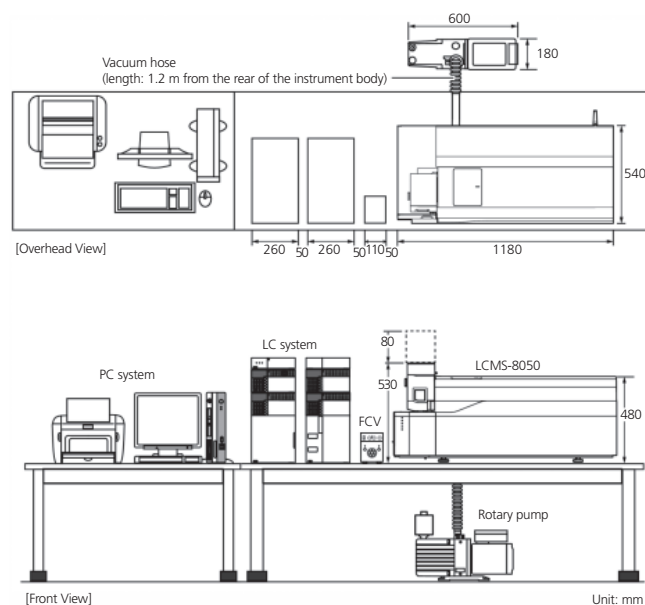

Installation Conditions

Temperature	18 to 28 °C
Humidity	20 to 70 % (Non-condensing)
Size	1180 mm (W) × 540 mm (D) × 610 mm (H)
Weight	140 kg
Power supply	MS unit: AC 230 V 15 A (50/60 Hz) Single-phase
Gas requirements	Nitrogen gas: Maximum 24.4 L/min, purity greater than 97 % Dry air : Maximum 20 L/min, oil/water-free Total nitrogen plus air: 25L/min maximum

The above are not standard installation specification. All LCMS-8050 instruments will be installed and tested in accordance with standard performance tests as detailed in the Shimadzu document ZEBN-5913, Shimadzu High-Performance Liquid Chromatograph Mass Spectrometer LCMS-8050 Installation Procedures.

This product is for research use only, not for use in diagnostic procedures.



Shimadzu Corporation
www.shimadzu.com/an/

Company names, product/service names and logos used in this publication are trademarks and trade names of Shimadzu Corporation or its affiliates, whether or not they are used with trademark symbol "TM" or "®". Third-party trademarks and trade names may be used in this publication to refer to either the entities or their products/services. Shimadzu disclaims any proprietary interest in trademarks and trade names other than its own.

Shimadzu does not assume any responsibility or liability for any damage, whether direct or indirect, relating to the use of this publication. The contents of this publication are subject to change without notice.